

mectron


medical technology




easyjet pro

Bicarbonate jet polisher

easyjet pro

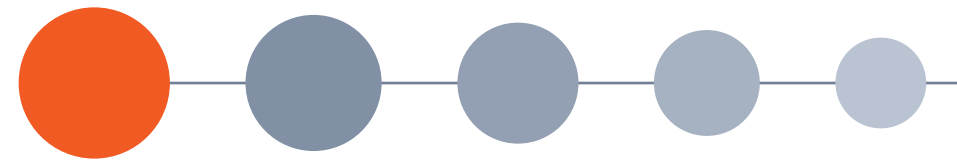

airpolishing for turbine coupling

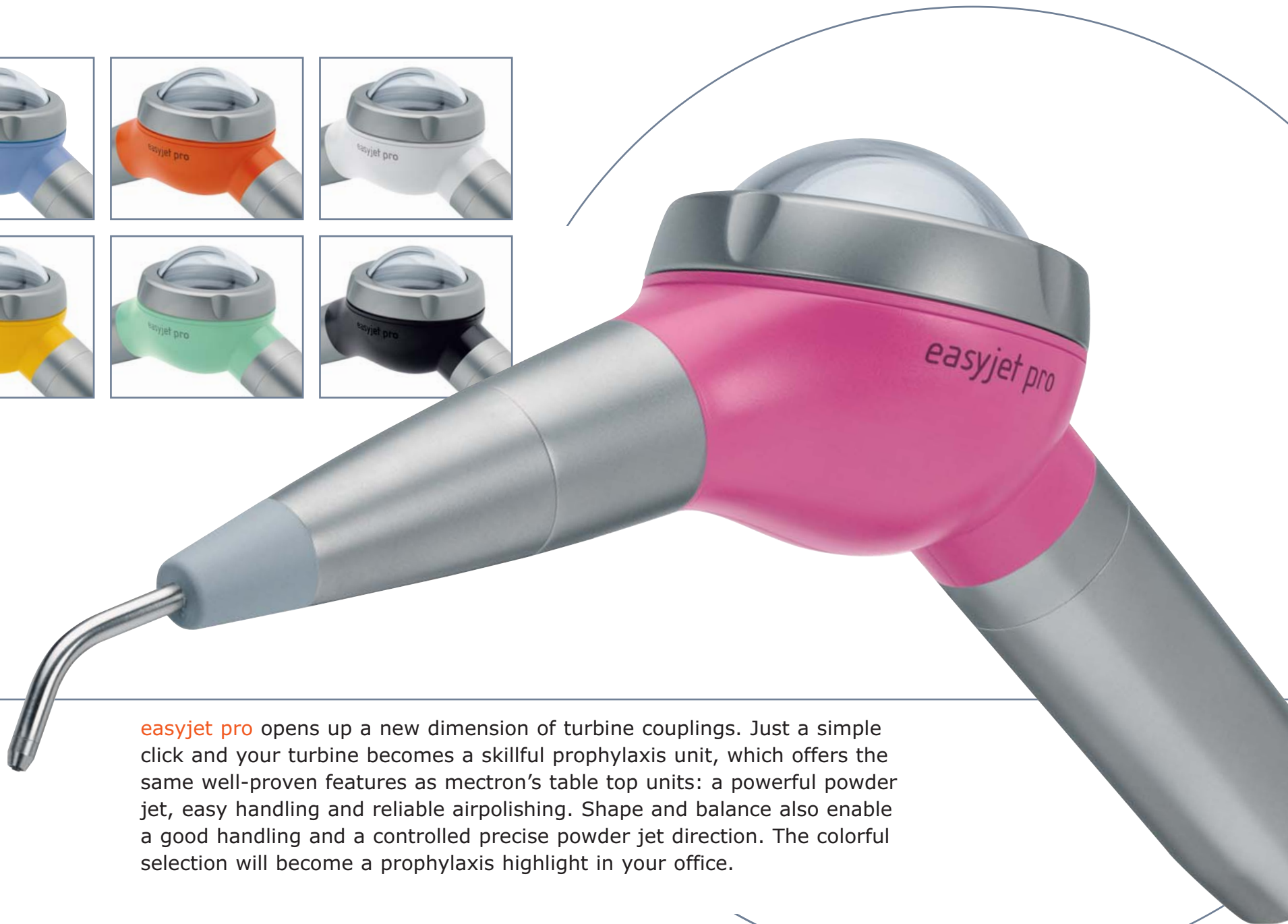

easyjet pro opens up a new dimension of turbine couplings. Just a simple click and your turbine becomes a skillful prophylaxis unit, which offers the same well-proven features as mectron's table top units: a powerful powder jet, easy handling and reliable airpolishing. Shape and balance also enable a good handling and a controlled precise powder jet direction. The colorful selection will become a prophylaxis highlight in your office.

Powder chamber

- powder capacity of 14 grams
- a new technique allows a constantly effective powder jet
- no clogging

Spray handpiece

- 120° angle
- just a click to connect
- 360° rotating head
- powder and water travel separately until they reach the nozzle
- the water surrounds the jet of powder and air and prevents powder dust
- fully sterilizable



Silicon protector

- protects teeth
- prevents any pollutant entry
- dismountable
- fully sterilizable



Powder chamber cap

- stainless steel ring-nut
- transparent dome for optimal control of powder level
- triangular thread for easy opening
- thanks to its specific shape there will be no powder deposits on the thread

Technical data

- incoming air pressure between 2.3 – 3.2 bar
- incoming water pressure between 0 – 3 bar




easyjet pro stands out due to its innovative technology, which allows high reliability. All parts are optimized for daily use.



the technology

reliable and ergonomic

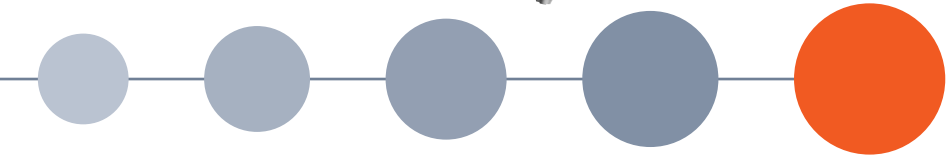

the turbine adapter

easyjet pro offers you a high degree of flexibility with its six different turbine coupling connectors. With just a click you can carry out the prophylaxis treatment. The connector can easily be changed by an authorized mectron service in case of a turbine coupling change.



Bien Air Unifix and Unifix L



KaVo multiflex and multiflex lux



Sirona R and B




W & H Roto Quick and Roto Quick lux



Midwest



Borden

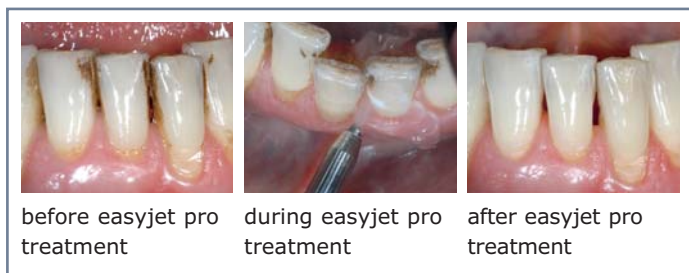

before easyjet pro treatment during easyjet pro treatment after easyjet pro treatment

Effect

Particles of sodium bicarbonate at a high rate of acceleration combined with warm water enable an extremely fine, regular and precise jet. The jet eliminates discolorations more effectively, gently and faster than when using rubber heads and polishing pasta, even in the interdental spaces.



Taste

The new taste is originated by a new lemon flavor mixed with a larger quantity of sweet saccharin.



Powder dosage

Each bottle has its dosage measuring cap, indicating the proper powder quantity for each mectron airpolisher.

Quality

Its optimized flowing and hydrophobic feature avoid clogging and allow a more regular powder jet with amazing treatment results.



the powder

mectron **prophylaxis powder** is delivered in convenient plastic bottles. Each bottle contains 250 grams of powder. The package includes two bottles of mectron **prophylaxis powder**.



MECTRON EASYJET PRO PLUG AND SPRAY!

 www.mectron.com or mectron@mectron.com

mectron s.p.a., via Loreto 15/A, 16042 Carasco (Ge), ITALIA, tel +39 0185 35361, fax +39 0185 351374



 **mectron**
medical technology

© Copyright mectron S.p.A., Carasco, Italy
All rights reserved. Texts, pictures and graphics of mectron brochures are protected by copyright and other protection laws. Without written approval of mectron S.p.A. the contents may not be copied, distributed, changed or made available to third parties for commercial purposes.